

S-Adenosyl-L-Methionine (SAM) ELISA Kit 3

Catalog No. Ik00203

Packaging size 48 tests

Detection range 7.5nM – 480nM

Target details

Methionine is an amino acid with sulfur methyl group, which is converted into S-adenosylmethionine (SAM) with ATP by methionine adenosyltransferase (EC 2.5.1.6). SAM is the sole methyl group donor for over 50 kinds of biologically active substances such as DNA, RNA, protein, phospholipids, hormones and neuro-transmitters, etc. After methyl group is transferred away from SAM under a methyltransferase, S-adenosylhomocysteine is formed, which is further metabolized into homocysteine (Hcy) after adenosine is removed. The homocysteine is metabolized to methionine through N5-methyl tetrahydrofolate methyltransferase (EC1.1.1.68) and coenzyme vitamin B12 catalytic pathway accepting methyl from N5-methyltetrahydrofolate. This is methionine cycle.

Methylation index is defined as the ratio of SAM and SAH. Methylation index is a better marker for methylation status and methylation capability.

SAM is used as medicines or nutritional supplements for treating depression, liver diseases, arthritis and joint pain, etc.

This immunoassay kit is used to measure the level of SAM in samples of soluble form.

Principle of the assay

This direct competitive ELISA (Enzyme Linked Immunosorbent assay) is designed to measure the level of S-adenosyl-L-methionine (SAM) in the sample. SAM conjugated with macromolecule is immobilized on the micro-titer plate. Standards and samples are pipetted into the wells, and then the HRP-conjugated antibody against SAM is added. The free SAM molecule in samples or standards competes with the immobilized SAM on the micro-titer plate surface for binding sites of the antibody. After discarding the mixed solution and washing each well, TMB substrate solution is added. The substrate solution turns blue under the effect of HRP (horseradish peroxidase), and changes into yellow once stop solution (acid) is added. The color develops in inverse proportion to the amount of SAM in the sample (or standards). The optical density of the remaining solution (OD450) is measured at 450nm using micro-plate spectrophotometer. The level of SAM in samples can be calculated through standard curve generated with standards.

Sample collection and storage

1. Sample collection must be carried out at 4°C. Remove precipitation by centrifugation when necessary, and test sample as soon as possible or store at -20°C or below. Avoid repeated freeze-thaw cycles.
2. Avoid Na₃N as it will inactivate HRP.

Procedure

1. Rewarm all reagents by putting them to room temperature and mix them well. Take out appropriate number of wells and put the remaining wells back into the original sealing bag. Seal and store it at -20°C.
2. Add 30ul sample diluent (blank well), standards and sample into each well.
3. Prepare HRP-antibody solution: Dilute HRP-antibody with HRP-antibody Diluent at 1:300, which should be used up within a week. Mix thoroughly and store in dark. Prepare the right amount that will be used in a week. Note: Always thoroughly centrifuge the HRP-antibody vial to recover all the 15-17ul HRP-antibody. It is recommended to collect all samples and measure them together in 1-3 experiments in a week by using 4.5ml of HRP-antibody Diluent to recover all HRP-antibody from the vial by washing it out several times.
4. Add 70ul HRP-conjugated antibody into each well except for the blank well.
5. Put the plate onto oscillator and shake for a while to mix the reagents, and seal the plate with micro-plate sealer. Incubate the plate at 37°C for 1 hour.
6. Peel the sealer carefully, and discard the remaining solution in the wells. Add at least 300ul wash solution in each well and maintain this state for 30 seconds, then remove wash solution. Repeat these steps 3 times to finish washing process. Or use auto-washer instead.
7. Add TMB Substrate and add 100ul blending substrate to each well. Shake gently and seal the plate. Incubate the plate at 37°C for 15 minutes without light.
8. Add 50ul Stop Solution to end reaction.
9. Measure absorbance of each well at 450 nm wavelength within 15 minutes. Set zero according to the blank well.

Data processing

Always clear the background by subtracting the absorbance at 450 nm (OD450) of blank well from that of the test well.

The binding rate of each well (standard or sample) is equal to A/AS_0 , with A being the average absorbance of the standard wells or the sample wells, AS_0 being the average absorbance of the S_0 standard wells.

Creation of standard curve: Construct a standard curve by plotting the binding rate of each standard on the y-axis against the logarithm of its concentration on x-axis. Use quadratic polynomial curve to fit the data ($r > 0.99$). The SAM level of the sample can be calculated by substituting its binding rate into the standard curve equation. Then multiply the extent of dilution if the sample has been diluted.

Note

1. Using reagents and samples without warm-up or ambient temperature less than 20°C may lead to reduced OD450 reading.
2. Extra drying of wells after washing may have negative effects on the results, such as poor standard curve and poor repeatability. To avoid it, perform the next step immediately after washing.
3. Mix solution well and wash completely, as these procedures will have influence on the assay.
4. Seal the plate with micro-plate sealer and avoid light when incubating the plate.
5. Duplicate wells for standards and samples are recommended, and quadratic polynomial is suggested for standard curve fitting ($R^2 > 0.99$). The detected concentration of quality control vial should be in the detection range.
6. The concentrated wash solution may form crystals. Warm it up to allow salts to dissolve completely before diluting.
7. The micro-plate sealers should be disposed after each use in order to avoid cross contaminations.
8. The HRP substrate (TMB) should be kept in dark.
9. The Stop Solution is diluted sulfuric acid. Avoid direct contact with skin and other things.

Storage and expiration

Storage condition: Except for the HRP substrate and stop solution that are stored at 2-8°C, all other ingredients and strips can be frozen stored.

Expiration: About 6 months from the date of shipment when storing it in refrigerator. If not planning to use it soon, users can put the assay plate, HRP-anti-SAM antibody, standards, HRP antibody and sample diluents in freezer for longer storage time or/and better results.

Components

Wash Solution (20×)	1×10ml
Assay Plate	6×8 wells (48-well)
HRP-antibody	1×15ul
HRP-antibody Diluent	1×4.5ml
Sample Diluent	1×0.5ml
TMB Substrate	1×5.5ml
Stop Solution	1×2.7-3ml
Microplate Sealers	1-2
S0 SAM Standard 0 nM	1×120ul
S1 SAM Standard 15 nM	1×120ul
S2 SAM Standard 30 nM	1×120ul
S3 SAM Standard 60 nM	1×120ul
S4 SAM Standard 120 nM	1×120ul
S5 SAM Standard 240 nM	1×120ul
S6 SAM Standard 480 nM	1×120ul
QC SAM Quality Control 160-240nM	1×120ul

Standard curve

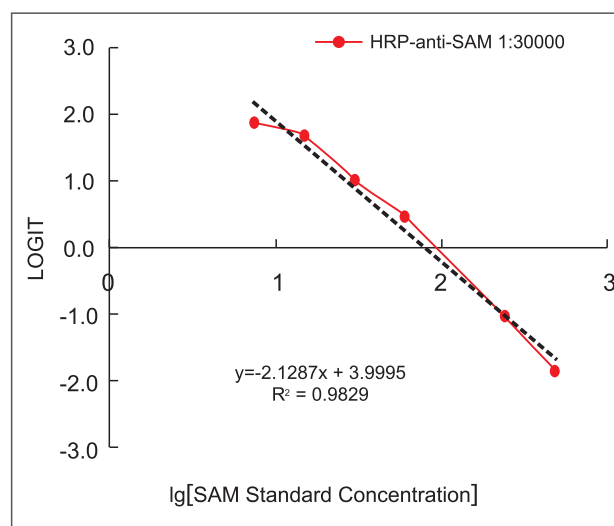


Figure 1 Standard curve HRP-anti-SAM antibody (clone 118-6) that was used at 1:30000 in a competitive ELISA. $LOGIT = \ln(A/AS_0/(1-A/AS_0))$. Sometimes adding 7.5nM to the standard impacts or slightly reduces the R of the standard curve. When the anticipated value of SAM from a sample is not less than 15nM, remove the 7.5nM point when plotting the standard curve. Further dilute the HRP-anti-SAM antibody to 1:35000 raises the R^2 to above 0.99 when the 7.5nM point is included.